

Instruction manual

Safe and efficient NiTi rotary system

Fulfilling the biological requirement for successful endodontics





Kit contents

BioRaCe kit

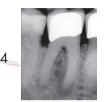


Table of contents

Concept	<u>4</u>
Minimal apical preparation sizes	<u>5</u>
Coding system	<u>8</u>
SafetyMemoDisc	<u>10</u>
Endo Stand	<u>11</u>
Basic Sequence	<u>12</u>
Extended Sequence	<u>17</u>
Obturation	<u>19</u>
Clinical cases	<u>20</u>
References	22







Pre-operative



Post-operative

Concept

Endodontics is the prevention or elimination of apical periodontitis

Since root canal infection is the cause of apical periodontitis, the biological aim of endodontic treatment is the prevention or elimination of root canal microbes.

Consistent success in endodontics requires high technical skill in order to achieve a biological aim. It is well established that in order to remove enough microbes from the root canal to ensure predictable success, the apical third of the canal must be instrumented to certain minimum sizes (see anatomical chart & references).


Most instrumentation systems require an additional step to achieve minimum sizes in the apical third of the canal. This results in additional files, time and expense for the practitioner.

The Bio*RaCe* sequence is unique, it has been especially designed to achieve the required apical sizes without the need for additional step and additional files. If used according to instructions, most canals can be effectively cleaned with 5 NiTi files. Thus with the use of the unique Bio*RaCe* system, the biologic aim of root canal treatment is achieved WITHOUT compromising efficiency.



Follow-up 1 year

Required minimal apical preparation sizes

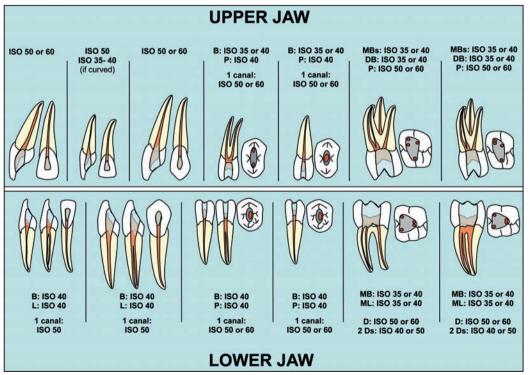

The attached chart describes the required minimal apical preparation sizes based on morphometric and anatomical studies of the root canal system (see references at the end).

Please, consult the anatomical chart before the treatment is initiated.






6

BioRaCe basic sequence

Note on the opposite figure that when the instruments are at full working length, the tip of BRO, BR1 and BR3 don't touch the canal walls (see contact zones).

Concept

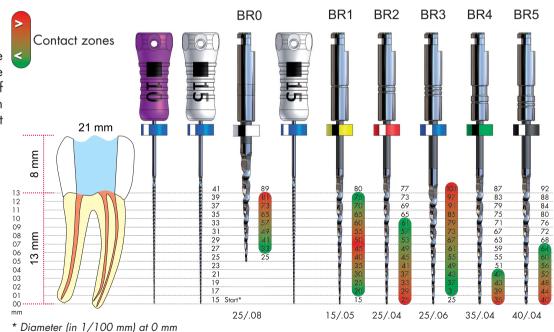

BioRaCe instruments present the same physical characteristics of RaCe instruments such as:



Bio*RaCe* differs from the well known *RaCe* instruments in regard to instruments sizes, tapers and sequence. The major goal of Bio*RaCe* is to achieve apical preparation sizes that are scientifically proven to effectively disinfect the canal (see references at the end).

Bio*RaCe* has been designed to clean the root canal efficiently and safely with few instruments.



BioRaCe should be run @ 500-600 rpm

Recommended torque : 1 Ncm

Rac



8

BioRaCe instruments are identified by a specific coding system on the handle and by specific colored rubber stoppers*

Example of marking on handles BRO to BR7 according to Basic sequence



Coding key (ck) Rubber Stopper (Assistant Colour Coding)




Coding system

Recommended torque value for all instruments:


1 Ncm

* The colored stoppers do not correspond to ISO color coding for sizes. However at the final apical sizes the colors are correct for ISO sizes (green #35 and black #40).

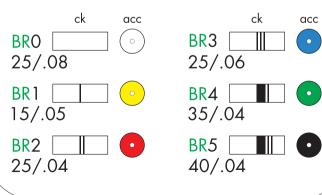

BioRaCe Basic Set

BioRaCe Extended Set instruments have an additional control device: SafetyMemoDisc*

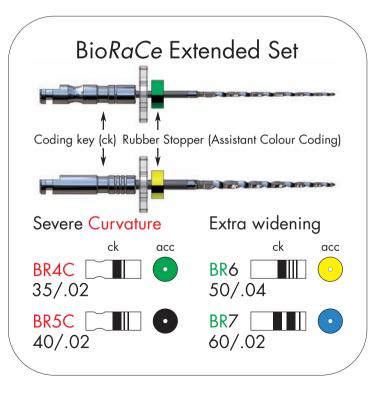



 Extra identification groove on
BR4C and BR5C instruments only (i.e. Severe Curvature)

Recommended torque value for all instruments:

1 Ncm



* With the SafetyMemoDisc (SMD) please see after








SafetyMemoDisc (SMD) device: a better control of risks for optimum safety

10

Increased safety:

- Discard instruments before they become a hazard; a new package is cheaper than spending time trying to remove broken pieces.
- Follow FKG recommendations on speed and torque to set your motor (see page 7).

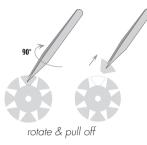

The SMD are mounted on the BR6/7/4C/5C instruments as a standard feature. You can control the fatigue of each instrument, according to the complexity of the canal.

Each *SMD* has eight petals. After each use, the practitioner will pull several petals off (see recommendations below). The remaining petals indicate the possibility of further use.

The SMD are sterilisable and therefore stay on the instruments. You keep the utilisation's information all along the file life.

How many times can BioRaCe BR6/7/4C/5C be used?

There is no straight answer – follow SafetyMemoDisc recommendations below.

SafetyMemoDisc

Example of canal complexity: Simple, radius ≥ 25mm Medium, radius 25 until 11 mm Difficult, radius ≤ 11 mm Remove: 2 petals 4 petals 6 petals



Endo Stand - dedicated instruments organizer



User-friendly, the BioRaCe Endo Stand is compact, strongly built to withstand all types of sterilizations, ensuring that the practitioner has the right tools on hand for carrying out the treatments. We recommend to use the Basic set instruments (BRO, BR1, BR2, BR3, BR4, BR5) in four cases, as a maximum. Therefore Bio*Race* Endo Stand has a special indicator of wear with 4 petals.

The practitioner will pull one petal off after each use. The remaining petals indicate the number of times the Basic set can still be used.

The information on the number of uses is recorded all along the Basic set life.



Indicator of wear, with 4 petals. Refill packs available separately



Note that there are 2 notches to help positioning the rubber stopper at the required length (suitable for lefthanded or right-handed persons)






Overview of the BioRaCe Basic sequence

12

The changes in sequence of sizes and tapers has allowed the required apical sizes to be achieved without increasing the number of instruments.

Basic

Sequence

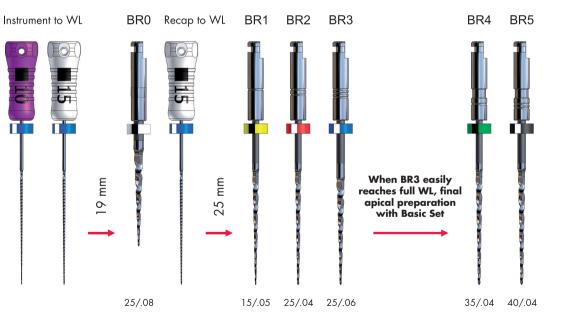


Pre-Operative Procedures

- Take parallel pre-operative radiograph(s).
- Place rubber dam (rubber dam may be placed after penetrating the roof of the pulp chamber, in case of a difficult access).
- Perform access.
- Localize canal with an endodontic probe.
- Remove coronal curvatures and establish straight line access to the canals orifices.
- Use disinfectant over the pulp chamber, tooth, and 1 cm of surrounding rubber dam.
- Establish working length with an Electronic Apex Locator (EAL) with SS k-files (e.g. 25 mm #08-15).





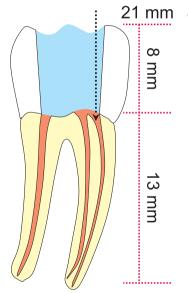







14

Manual Instrumentation Phase

SSt files to Working Length (hereafter WL) # 08, 10, 15

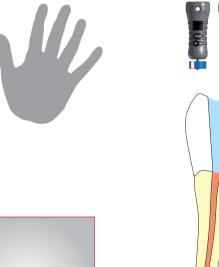

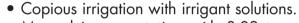

- Manual instrumentation with 0.02 taper SSt files from # 08 to #15 to full WL.
- Irrigation.

8 mm


13 mm

......

• If needed verify radiographically the WL with a SSt file #15.

Instrumentation phase

Rotary Instrumentation Phase: Access by using BRO instrument

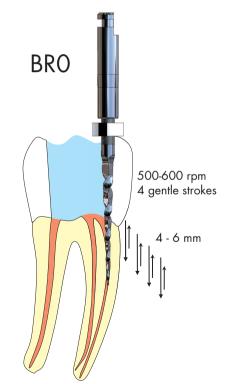

- Do not start this phase until a K-file #15 comfortably reaches the WL.
- Adjust motor to 500-600 rpm and 1 Ncm.
- Fill the canals and pulp chamber with irrigant.
- BRO "only" 4 gentle strokes clean the flutes.
- Repeat until approximately 4-6mm of coronal part of the canal has been prepared.









Rotary Instrumentation Phase: Reach Working Length (WL) with BR1 to BR3



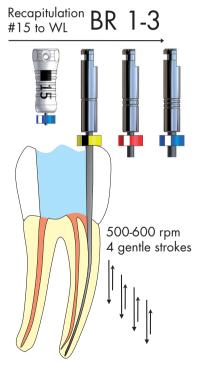

- After use of BRO, repeat irrigation.
- Recapitulate to full WL with a SSt file #15.
- Fill the canal and pulp chamber with irrigant.
- Use BR1 with 4 gentle strokes. If this instrument does not reach the WL, clean the instrument and repeat until the WL is achieved (If necessary, reconfirm the WL with an Electronic Apex Locator).
- Use BR2 and BR3 as described for BR1.
- DO NOT use BR3 to full WL on canals with severe apical curvatures.
- Irrigate copiously between instruments.

Rotary Instrumentation Phase: Final apical preparation with BR4 to BR7

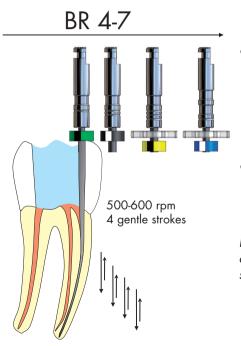

- In most cases, the final apical preparation is achieved with instruments BR4 and BR5. Depending on the root canal anatomy (see anatomical chart), two additional instruments BR6 and BR7 can be used for larger canals.
- The same principle as explained for BR1-3 should be used for the apical preparation.

NB. Copious irrigation at all times and cleaning of the files after 4 gentle strokes is essential for safe and efficient use of these instruments.







Cases with severe apical curvatures Specific instruments: BR4C and BR5C



Instrumentation phase


BR 4C-5C 500-600 rpm 4 gentle strokes

- For severe apical curvatures, instruments BR4C and BR5C should be used to prepare the apical canal.
- If the instrument does not reach the WL with 4 gentle strokes, DO NOT FORCE the instrument. Irrigate the canals and repeat.
- For complicated curvatures it is recommended to use additional FKG instruments (e.g. S-Apex inverted taper instruments).





Obturation phase

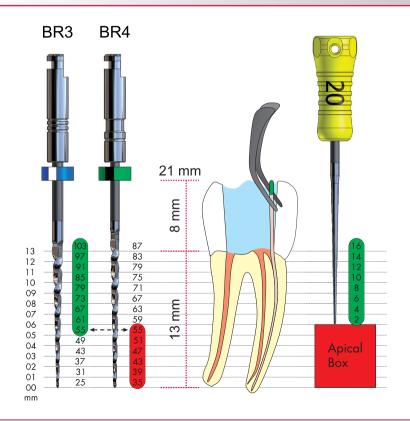

- Place obturation point corresponding to the final apical sizes, in this case BR4 35/0.04.
- Complete with preferred obturation material and technique.
- For lateral condensation technique a NiTi finger spreader ISO 20/0.04 is suggested.
- In cases where a warm vertical condensation technique is planned to be used, a fine or a fine-medium heat plugger should be used when the final apical preparation has been achieved with BR4 or BR5.





20

Clinical cases: Moderate curvatures

Tooth 25 Dx: Symptomatic Pulpitis Tx: Pulpectomy

More world wide clinical cases on www.biorace.ch Treatment Details: MB #35/0.04 DB #35/0.04 D #50/0.04

Tooth 46 Dx: Asymptomatic apical periodontitis Tx: Of non vital tooth

Treatment Details: MB #35/0.04 ML #35/0.04 DB #50/0.04 DB #50/0.04







Clinical cases: Severe curvatures

Tooth 26 Dx: Symptomatic Pulpitis Tx: Pulpectomy

Treatment Details: MB1 and 2: #35/0.04 DB #40/0.04 P #60/0.02

Tooth 27 Dx: Asymptomatic apical periodontitis Tx: Of non vital tooth

Treatment Details: MB1 and 2: #35/0.04 DB #40/0.04 P #50/0.04





More world wide clinical cases on www.biorace.ch





References

22

References

- Bartha T, Kalwitzki M, Löst C, Weiger R, Extended apical enlargement with hand files versus NiTi files. Part II. Oral Surg Oral Med Oral Pathol Oral Radiol Endod 2006; 102: 692-7.
- 2) Baugh D, Wallace J. The role of apical instrumentation in root canal treatment: a review of the literature. J Endod 2005; 31: 333-340.
- 3) Byström A, Happonen R, Sjögren U, Sundqvist G. Healing of periapical lesions of pulpless teeth after endodontic treatment with controlled asepsis. Endod Dent Traumatol 1987; 3: 58–63.
- 2) Card SJ, Sigurdsson A, Orstavik D, Trope M. The effectiveness of increased apical enlargement in reducing intracanal bacteria. J Endod 2002; 28: 779–783.
- 3) Dalton BC, Ørstavik D, Phillips C, Pettiette M, Trope M. Bacterial reduction with nickel-titanium rotary instrumentation. J Endod 1998; 24: 763–7.
- 4) Gomes BP, Souza SF, Ferraz CC, Teixeira FB, Zaia AA, Valdrighi L, Souza-Filho FJ. Effectiveness of 2% chlorhexidine gel and calcium hydroxide against Enterococcus faecalis in bovine root dentine in vitro. Int Endod J 2003; 36: 267–275.
- 5) Kerekes K, Tronstad L. Morphometric observations on root canals of human anterior teeth. J Endod 1977; 3: 24–29.
- 6) Kerekes K, Tronstad L. Morphometric observations on root canals of human premolars. J Endod 1977; 3: 74–79.
- 7) Kerekes K, Tronstad L. Morphometric observations on the root canals of human molars. J Endod 1977; 3: 114–118.
- 8) Kerekes K, Tronstad L. Long-term results of endodontic treatment performed with a standardized technique. J Endod 1979; 5: 83–90.
- 9) McGurkin-Smith R, Trope M, Caplan D, Sigurdsson A. Reduction of intracanal bacteria using GT rotary instrumentation, 5.25% NaOCl, EDTA, and Ca(OH)2. J Endod 2005; 31: 359–63.
- 10) Mickel AK, Chogle S, Liddle J, Huffaker K, Jones JJ. The role of apical determination and enlargement in the reduction of intracanal bacteria. J Endod 2007; 33: 1:21-23

References

- Safavi KE, Nichols FC. Effect of calcium hydroxide on bacterial lipopolysaccharide. J Endod 1993; 9: 76–78
- 12) Shuping GB, Ørstavik D, Sigurdsson A, Trope M.Reduction of intracanal bacteria using nickeltitanium rotary instrumentation and various medications. J Endod 2000; 26: 751–755.
- 13) Siqueira, J.F., Jr & de Uzeda, M. Disinfection by calcium hydroxide pastes of dentinal tubules infected with two obligate and one facultative anaerobic bacteria. Journal of Endodontics 1996; 22, 674-676.
- 14) Sjögren U, Figdor D, Persson S, Sundqvist G. Influence of infection at the time of root filling on the outcome of endodontic treatment of teeth with apical periodontitis. Int Endod J 1997; 30: 297–306.
- 15) Sjögren U, Figdor D, Spångberg L, Sundqvist G. The antimicrobial effect of calcium hydroxide as a short-term intracanal dressing. Int Endod J 1991; 24: 119–125.
- 16) Spängberg L, Rutberg M, Rydinge E. Biologic effects of endodontic antimicrobial agents. J Endod 1979; 5: 166 –75.
- 17) Teixeira FB, Levin LG, Trope M Investigation of pH at different dentinal sites after placement of calcium hydroxide dressing by two methods. Oral Surgery, Oral Medicine, Oral Pathology, Oral Radiology, and Endodontology 2005; 99: 511–6.
- 18) Trope M, Debelian G. Endodontic treatment of apical periodontitis. in: Ørstavik D, Pitt Ford T; Essential Endodontology; 2nd; Blackwell, Munksgaard 2007.
- 19) Trope M, Debelian G. Endodontics manual for the general dentists. Quintessence publishing, UK 2005. also translated to Polish, Russian and Turkish.
- 20) Weiger R, Bartha T, Kalwitzki M, Löst C. A clinical method to determine the optimal apical preparation size. Part I. Oral Surg Oral Med Oral Pathol Oral Radiol Endod 2006; 102: 686-91





Safe and efficient NiTi rotary system

Your FKG's partner:

For more information, please view supplied Educational CD or visit www.biorace.ch



FKG REF. - 99.7AA.00.02D.AN/Rev. 4 - 08/2012