

# Scout-RaCe


## Mechanical scouting sequence

### Description

*Scout-RaCe* enables a faster and better shaping of glide path than manual files. *Scout-RaCe* 10/.02 can be used directly after a manual instrument ISO 08 up to working length, followed by *Scout-RaCe* 15/.02 and *Scout-RaCe* 20/.02. Thanks to its exceptional flexibility, *Scout-RaCe* perfectly follows the canal anatomy.

*Scout-RaCe* offers the following features:

- Extreme flexibility thanks to the .02 taper
- Exclusive rounded Safety Tip
- Unique anti-screwing design
- Sharp cutting edges for efficiency
- Electro-chemical polishing: better resistance to torsion & fatigue
- *SafetyMemoDisc*: for mastering metal fatigue & controlling number of uses

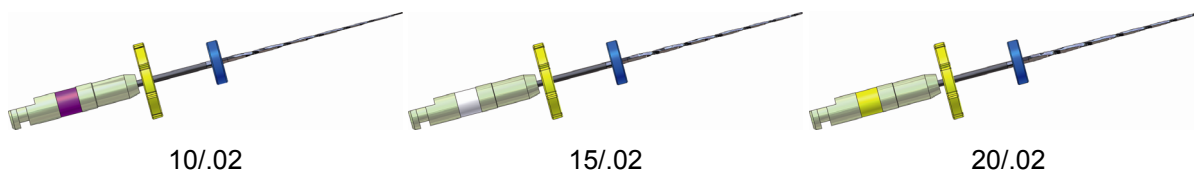

### Availability

- Assortment (2 pces of each)
- Individual sizes (5 pces/box)