



# BC UNIVERSAL SEALER & RRM

# BC UNIVERSAL SEALER




**Easy-to-use premixed calcium aluminosilicate paste**

**Indication:**  
Permanent obturation of the root canal


## Biocompatible

- High purity calcium aluminosilicates
- Resin-free
- Eugenol-free


## High pH > 12

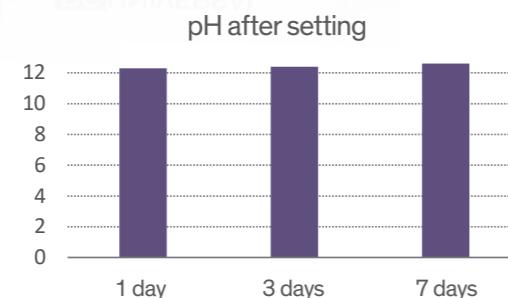

## Radiopaque

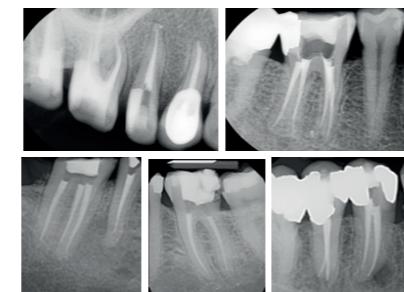


## Hydroxyapatite generation

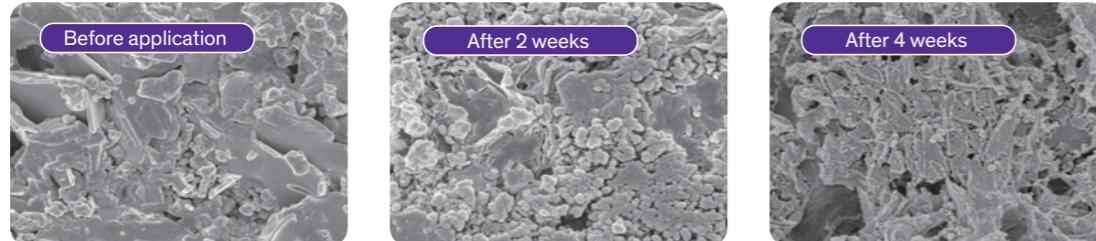

## Excellent Root Canal Sealing

- Excellent bonding to gutta-percha and dentin
- Zero shrinkage
- Optimal flowability

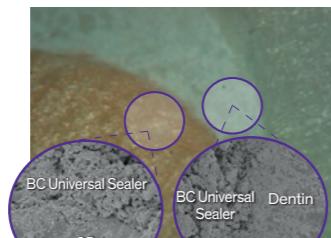

Optical microscope and SEM image of GP/BC Universal Sealer/Dentin interface : No gap at interface



Optical microscope image of microleakage test (dye penetration test):  
No dye penetration

## Product reference

## Content

## Reference Number



BC Universal Sealer

2g Syringe + 20x Disposable tips + Syringe Holder

08.650.00.013.XX


# BC UNIVERSAL RRM



**Premixed root canal filling material  
Calcium aluminosilicate paste  
(putty type)**


**Indications:** Pulp capping, Repair of root perforation, Repair of root resorption, Root-end filling, Apexification

## Highly Biocompatible

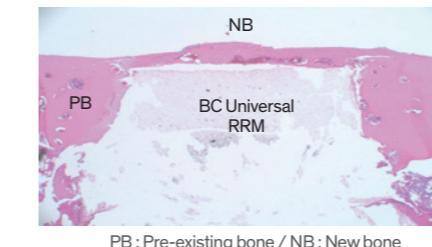

PB : Pre-existing bone / NB : New bone  
BC Universal RRM does not create inflammatory response


## High pH ≥ 11

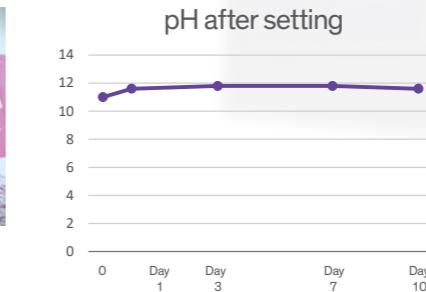

- Excellent compressive strength (around 100 MPa)
- Superior dislodgement resistance due to excellent bonding to dentin
- High stability due to good hydration reaction
- Low solubility
- Excellent color stability
- Appropriate setting time for application
- Radiopaque

## High mineralization ability

- BC Universal RRM forms extensive apatite crystals on the surface
- BC Universal RRM provokes better cell adhesion and cell spreading vs. white MTA

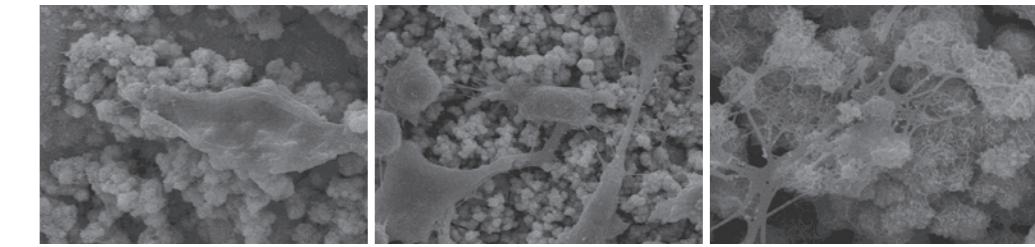

White MTA (X2,000)

BC Universal RRM (X2,000)

BC Universal RRM (X10,000)

## Product reference

## Content

## Reference Number



BC Universal RRM

0.25g Capsule x10

08.650.00.014.XX

## Product reference

## Content

## Reference Number



BC RRM

0.25g Capsule x10

08.650.00.014.XX



# Empower Endo

For 90 years, we have striven to deliver Swiss quality, precision and peace of mind to endo specialists. Our unique approach is driven by the quest to provide the best dental care and our conviction that

the search for perfection is never over. Like you, we seek to advance the frontiers of the possible by evolving through research, refinement and innovation. Ready to change the game?