



FKG
swiss endo



NiTi rotary files
for mechanised
glide path preparation

EN

ScoutRace

Sequence designed for canals with severe curvature



► The FKG Race range

An internationally-recognised reference, Race (Reamer with Alternating Cutting Edges) is the most complete selection of NiTi rotary files for root canal treatment.

ScoutRace

The range designed for canals with severe curvature

ScoutRace is a sequence of three instruments with .02 taper and with ISO sizes of 10, 15 and 20. They are designed to mechanically prepare root canals with severe curvature, or «S» type canals. The files operate at a speed of 600-800 rotations per minute, and are used after determining the working length. The ScoutRace range enables the glide path shaping faster and safer than

with manual instruments. The anatomy of the median and apical thirds are well respected, to experience no transportation and nor ledges. The files adapt perfectly to the canal's shape, thanks to their .02 taper and remarkable flexibility. The ScoutRace can be used without pressure up to working length. Root canal shaping can then be carried out with dedicated FKG instruments.

Illustration of a case with a severe apical curvature

Courtesy of Dr. Gilberto Debelian (Norway)



Pre-operation
Diagnosis: pulpitis

Difficulty:
Double curve

Post-operation
Preparation: 35/.04

FKG Race: Specifications


Several exclusive technological advances lie behind Race range's efficacy and worldwide renown:

- Exclusive rounded safety tip for precise guiding
- Alternating cutting edges to eliminate screwing-in effect
- Sharp edges for best cutting efficiency
- Electro-chemical polishing for better resistance to torsion and fatigue
- SafetyMemoDisc (SMD) to master metal fatigue

ScoutRace availability

Sterile blister 25 mm
2 × 1 ScoutRace Sequence
6 × Individual files

Non-sterile blister 21/25/31 mm
2 × 1 ScoutRace Sequence
6 × Individual files

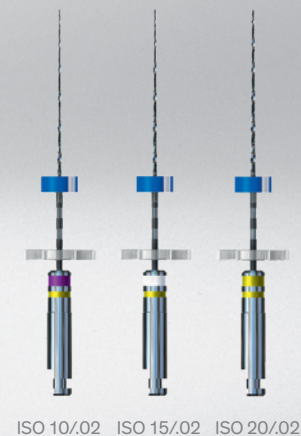



FKG Dentaire Sàrl

Founded in Switzerland in 1931, FKG Dentaire Sàrl gained a new momentum in 1994, the year Jean-Claude Rouiller took over the reins of the company.

He propelled FKG to the forefront in the development, manufacturing and distribution of dental products destined for general practitioners, endodontists and laboratories.

The FKG strategy is centered on innovative high-precision products and the creation of machines designed specifically for the dental field. Its aim is to offer solutions that meet the most demanding needs of end users.

In 2011 the son of Jean-Claude Rouiller, Thierry, succeeded to the head of the company. Through his

incentive, the network of distributors has expanded significantly and allowed FKG to make its products available in over 100 countries worldwide.

In 2012, the Swiss Venture Club rewarded FKG for its dynamism, high product quality, and its continuing innovation.

Equipped with a clean room since 2013, FKG is now developing a range of sterile products.

In 2013 and 2014 the company unveiled state-of-the-art training centers in La Chaux-de-Fonds, Dubai, and Oslo.

FKG Dentaire is certified according to international norms and regulations.



FKG Dentaire Sàrl, Le Crêt-du-Loche 4, 2322 Le Crêt-du-Loche, Switzerland
T +41 32 924 22 44, info@fkg.ch, www.fkg.ch