

D-RaCe


New desobturation instruments

"Desobturation" prior to "retreatment": because access must be made, the entrance of the canal "desobtured" to reach inside the canal itself so that it could be properly retreated. Retreatment must be achieved because the first RCT failed or needs improvement.

Description


D-RaCe will desobturate most canals, whatever the filling used – gutta-percha, carriers, paste, resin-based materials. D-RaCe removes a large part of the ancient filling material, and without solvent. With only 2 instruments time is saved and safety of treatment is increased.

Sequence at a glance

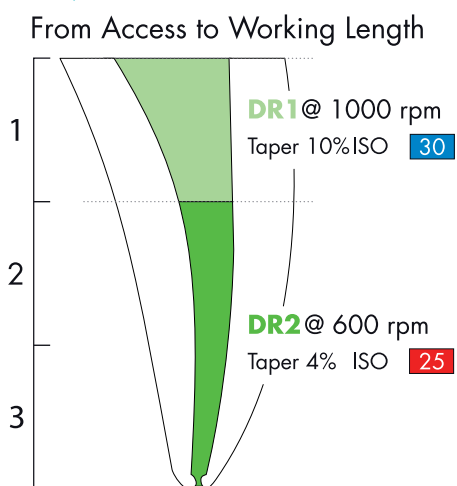

DR1 - Access

- >>> ISO 030/0.10 - L.15/8 mm
- >>> active tip

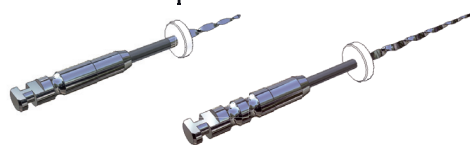

DR2 - Full path

- >>> ISO 025/0.04 - L.25/16 mm
- >>> safety tip

Availability

- Assortment (2 pces of each)
- Individual sizes (6 pces/box)

Website

www.d-race.ch