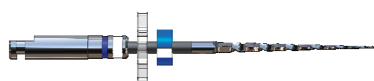


**FKG Dentaire introduces a new quick, effective and safe sequence: *iRaCe*****Description**


Thanks to their exclusive features, only 3 *iRaCe* rotary NiTi files are needed to treat most cases (straight, slightly curved and/or large). *iRaCe* sequence is easy to learn and to use, meaning considerable time saving.

The *iRaCe* rotary NiTi files offer the following features:


- *Exclusive rounded Safety Tip*: Perfect guidance
- *Alternating cutting edges*: No screwing-in effect, no pulling in
- *Sharp cutting edges*: They cut efficiently = saves time
- *Electrochemical polishing*: Better resistance to torsion and cyclic fatigue
- *SafetyMemoDisc (SMD)*: To master fatigue and number of uses
- *New shank*: Easy identification of ISO sizes (large ring, yellow: 2%, red: 4%, blue: 6%)

*iRaCe*

15/06

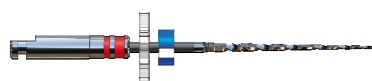


25/04

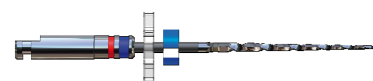


30/04

Two additional instruments have been developed to allow the treatment of difficult cases (severely curved, narrowed and/or calcified): *iRaCe Plus*

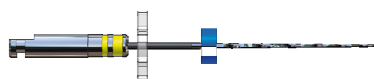*iRaCe Plus*

20/02

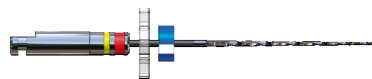

25/02

**Availability**

- *iRaCe* Sequence (3 instruments/blister)
- *iRaCe Plus* (2 instruments of each/blister)
- Refill (6 pcs/blister)
- Length 21 / 25 / (31 mm, longer delivery time)

**Website**[www.irace.ch](http://www.irace.ch)